

Desktop Clean Hood PVC 탁상용 클린 후드

배송비 별도





Cat. No.	Model	Size	EA Price ₩
C10-91-201	LM-0700P	w700×d720×h1000 mm	3,343,400
C10-91-203	LM-0900P	w900×d720×h1000 mm	3,893,800

- 우기미량 분석용 시험장비의 준비작업 등 청정공간을 유지하여 약품 증기 유출을 차단하는 장비
- 청정 공간에서 산 작업이 필요한 경우에 사용
- 공기가 상부에서 필터를 거쳐 수직 방향으로 하강하여 후면으로 배기 되는 방식
- 내부가 비금속 재질로 구성되었기 때문에 부식에 의한 오염이 없음
- 균일한 Lamina Flow 를 형성하여 신뢰성을 높임
- Class 100 HEPA Filter를 사용하여 99.97%의 집진율
- FL Lamp 및 UV Lamp 적용
- Filter 교환시기를 경보음과 LCD 화면에 표시
- 차압계로 필터의 상태를 확인 할 수 있음
- 10단계 풍속 조절 가능
- 2단 접이식 홀딩 도어
- Timer 기능 (99 hour 59 min)
- **PL Insurance (제조물 배상 책임 보험)**
- Validation Service (IQ, OQ, PQ)
- LK Labkorea

Cat. No.		C10-91-201	C10-91-203
Model		LM-0700P	LM-0900P
Size	Internal	w680×d570×h590 mm	w880×d570×h590 mm
	External	w700×d720×h1000 mm	w900×d720×h1000 mm
	Exhausted Duct	L150, d200 mm	
Filter	HEPA	HEPA filter 0.3 um 99.97% Efficiency Filter	
	Pre-filter	Efficiency Nylon Filter	
Material	Upper Case	Powder Coated Steel	
	Working Space	PVC	
Wind Velocity		0.3 ~ 0.5 m/sec, 10 step adjustable	
Lamp		FL lamp 1ea, UV lamp 1ea	
Control System		LCD Digital Controller	
Power		1 Phase, 220 Vac, 60 Hz	

Optional Accessory

Cat. No.	Model	Description	EA Price ₩
C10-91-903	LM-FL10	FL lamp 10 W	8,300
C10-91-913	LM-UV10	UV lamp 10 W	30,000
C10-91-940	LM-CHF700	Filter for LM-0700P	217,200
C10-91-941	LM-CHF900	Filter for LM-0900P	312,800

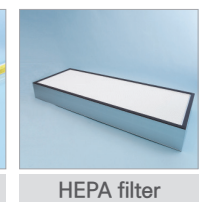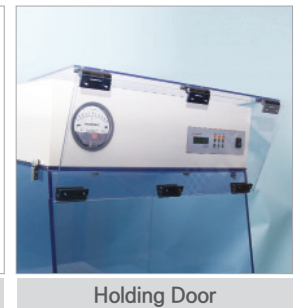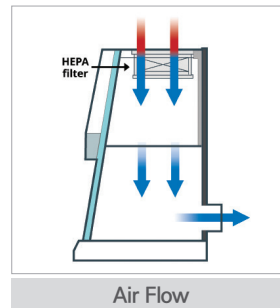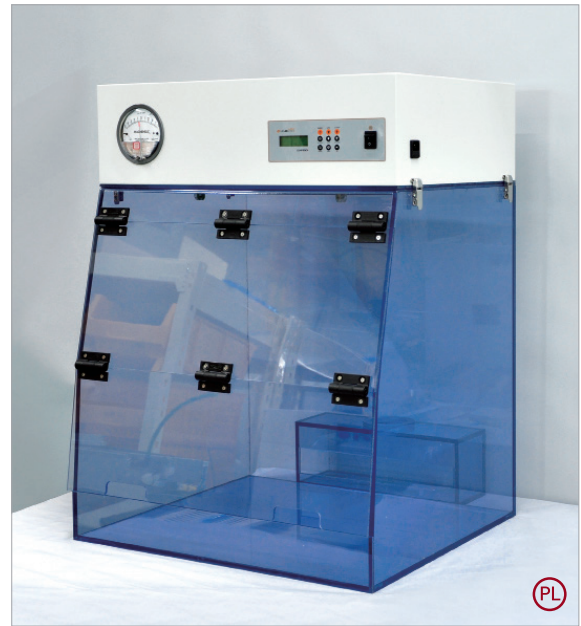

- 필터식 Fume Hood는 626 Page 참조
- 필터식 Cabinet은 528 Page 참조

