

## Cold Trap Bath 저온 동결 트랩

배송비 별도

Cat. No.	Model	Capacity	Temp. Range	EA Price ₩
B04-01-230	CTB40	4.2 L	Down to -40°C	2,682,800
B04-01-240	CTB80		Down to -70°C	4,813,900

- 진공 라인에서 발생하는 증기 및 유해성 물질 응축용
- 온도 조절 불가능
- 설치가 용이한 2개의 Glass Trap 사용 (기본 포함)
- 컴팩트 저소음형 동결 트랩
- 착탈식 전면 커버 구조, 냉동기 유지 관리 용이, 효율 증대
- Bath 내부 용액을 편리하게 배출 시킬 수 있는 Drain 장착
- Graphic LCD Display, 간편하고 직관적인 Controller 조작
- 온도 Calibration 가능
- Automatical Fault Sensing System 기기 스스로 이상 상황을 감시 (특허 등록)
- Validation Service (IQ, OQ, PQ)
- **PL Insurance (제조물 배상 책임 보험)**
- LK Labkorea

### Optional Accessory

Cat. No.	Model	Description	EA Price ₩
B04-01-A30	GCT01	Spare Glass Trap	52,000
B04-01-A40	GCT02	Spare Rubber Gasket	65,800

### Validation Service

Cat. No.	Model	Description	EA Price ₩
B04-91-995	VA-CTB100	Validation for Cold Trap	

- 당사 기본 양식에 따른 IQ/OQ Validation Service 비용이며, 지역이나 고객 요청사항에 따라 비용이 변경될 수 있음




Drain



Bath Cover



Glass Trap

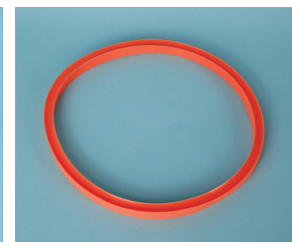

Rubber Gasket

Cat. No.	B04-01-230	B04-01-240	
Model	CTB40	CTB80	
Capacity	4.2 L		
Controller	Digital Control		
	Resolution 0.1°C		
	Timer Run start / Wait start (99 hr 59 min)		
Temperature	Range Down to -40°C	Down to -70°C	
	Accuracy ± 0.3°C		
Size	Interior Ø210×h220 mm		
	Exterior w340×d440×h640 mm	w450×d610×h750 mm	
Material	Interior Stainless Steel 304		
	Exterior Powder Coated Steel		
	Cover PMMA		
	Gasket Silicon Rubber		
Electric Supply	Power 1 Phase, 220 Vac, 60 Hz		
	Max Consumption 400 W (1.9 A)	800 W (3.7 A)	
	Power Line Standard Plug		
Other	Drain 1/4" Hose Nipple Male		
	Caster Leveling Caster		
Safety	Over Current Protection, Automatical Fault Sensing System		

Vacuum Oven은 744 ~ 748 Page 참조

