

Digital Sieve Shaker, Octagon® 디지털 체 진동기, up to Ø203mm Sieve

Cat. No.	Model	사용 sieve	EA Price ₩
S20-31-600	Octagon D200	Ø75, 100, 200, 203 mm Sieve용	6,609,600

- 고성능의 디지털 시브 셰이커 • 최대 직경 8" Sieve까지 사용 가능한 Shaker
- 20 um 이하의 Sieving을 위해 Wet Sieving Kit을 사용 가능
- Quiet, Electromagnetic Vibration Action
- 3 D Performance : Vertical Vibration & Rotating Action
- Shaker의 수직운동은 Sieve의 구멍이 막히는 것을 막아줌
- Start / Stop 시간을 설정 가능
- There are no rotating parts consequently
it is quiet in operation and maintenance free
- Anti-vibration Feet • Test Sieve는 별도 구매
- Complies with the Requirements of AASHTO T27 • Endecotts / Octagon®

Cat. No.	S20-31-600
Model	Octagon D200
Vibration Speed	3000 rpm
Approx. Amplitude	3.0 mm
최대 거취 수량	8개 (Ø200 ~ 203mm), 16개 (Ø76 ~ 100mm)
Timer Mode	Intermittent
Intermittent Process Timer	0 to 99.9 minutes
Intermittent Cycle Timer	Off time: 1 to 99 Seconds on time: 1 to 99 Seconds
Display	Digital 0.3 LED
Memory	Non Volatile, One set of data Stored
Size	w418 × d435 × h232 mm
Weight	35 kg
Power	230 V, 50/60 Hz
Class	I (Earthed)


Wet Sieving Conversion Kit, 습식 작업용 키트

Cat. No.	Model	Description	EA Price ₩
S19-251-920	RW-200	Wet Sieving Conversion Kit for 200mm	1,910,800
S19-251-930	RW-203	Wet Sieving Conversion Kit for 203mm	1,910,800

- 구성 : Top Clamping Plate with Perspex Cover, Spray Rose, Water Tight Seal, Stainless Steel, Receiver with Drainage Spout



S20-31-600

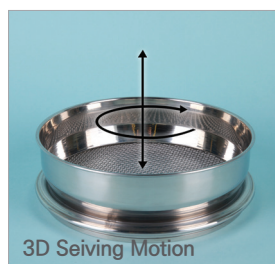

3D Seiving Motion



Receiver

Compact Sieve Shaker, Minor® 경제형 체진동기, up to Ø203mm

Cat. No.	Model	사용 sieve	EA Price ₩
S20-31-250	Minor M200	Ø75, 100, 200, 203 mm Sieve용	3,101,200

- 경제형으로 간단하게 사용할 수 있는 체진동기
- Compact 형으로 이동이 간단함 (17 kg)
- 최대 직경 8" Sieve까지 사용 가능한 Shaker
- 75 mm, 100 mm, 200 mm와 203 mm (8") Sieve용 Shaker
- 진동 고정식 Sieve Shaker • Quiet, Electromagnetic Vibration Action
- 3 D Performance : Vertical Vibration & Rotating Action
- Shaker의 수직운동은 Sieve의 구멍이 막히는 것을 막아줌
- There are no rotating parts consequently
it is quiet in operation and maintenance free
- Wet Sieving Operation용으로는 부적합 • Anti-vibration Feet
- Test Sieve는 별도 구매
- Complies with the Requirements of AASHTO T27 • Endecotts / Minor®

Cat. No.	S20-31-250
Model	Minor M200
Vibration Speed	3000 per min
최대 거취 수량	8개 (Ø200 ~ 203mm Sieve), 12개 (Ø100mm Sieve)
Timer	Continuous / 0 ~ 60 min
Size	Ø262mm(하부), h126mm(+Clamp Plate 30mm)
Weight	17 kg
Power	230 V, 50/60 Hz
Class	I (earthed)



S20-31-250